



medical **ECONET**
GERMANY

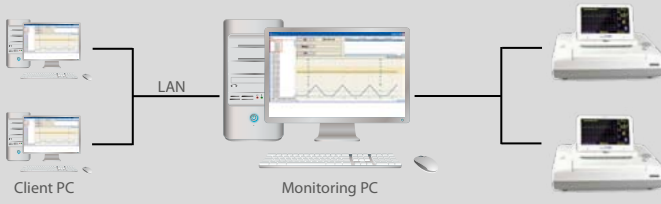


ECOtwin LCD

Modernized user-friendly antepartum fetal monitor

- 7" Color Display (800 x 480), rotatable by 270°
- 9 - crystal US-probes with 0,985 MHz frequency
- Waterproof IXP8 (US & UC)
- Integrated Li-Ion battery
- Two USB ports for data storage
- LAN port for data transmission
- Short Term Variability (STV) Analysis
- CCV (Cross-Channel Verification) Function
- 450 hours patient data saving
- Incl. probes, event marker, CTG-belts, ultrasound gel and paper
- Multiple language support
- Free Central Software

Specifications

Ultrasound
Multicrystal transducer with 9 elements
Pulsed doppler
Ultrasound frequency: 0.985 MHz
Intensity: <math><10\text{mW} / \text{cm}^2</math>
FHR range: 30 ~ 240 bpm
FHR accuracy: $\pm 2\%$ of range
Waterproof IPX 8
UC (Uterine Contraction)
External type
Reference (Zero) Control
Measurement range: 0 ~ 99
Waterproof IPX 8
FM (Fetal Movement)
Detection method: Ultrasound
Dual FM
Display
7" wide TFT color LCD
Heart rhythm
Printing state
Alarm on / off state LED
CTG Analysis
Baseline, OSC. Amp., OSC Num.
Acceleration, Deceleration (early, late, variable)
Movement
Short-Term VAR (min, max, mean)
FREE Central Monitoring System
Up to 8 fetal monitors
Recording and playback of Fetal Heart Rate
Multi-language support
Networking via LAN


Printer
Thermal Z-fold type (width: 150 mm)
Record speed: 1, 2, 3 cm / min
Record contrast: mode 1, 2, 3
Auto record period: Off, 10, 20, 30, 40, 50, 60 min
Set-up
Alarms upper / lower value
Time & Date
Print Contrast 1, 2, 3
Print speed
FHR volume level
CCV (Cross-Channel Verification)
Function
Mark function
UC reference
Zoom in particular response
Signal Interface
2x USB, LAN, RS-232
Environmental Specifications
Temperature Range
Operating: 10 to 45°C
Storage: -20 to 60 °C
Relative Humidity Range
Operating: 5 ~ 85 %
Storage: 20 ~ 95 %
Power Specifications
Power Requirements: 100 - 240 VAC, 50 / 60 Hz
Output: DC + 18V 2.8 A
Battery: Rechargeable, 14.8 V Li-Ion
Optional Accessories
Mobile trolley



CE 0123



Im Erlengrund 20, 46149 Oberhausen - Germany

t: +49 (0) 208 377 890-0 info@medical-econet.com
 f: +49 (0) 208 377 890-55 www.medical-econet.com
 facebook.com/medicaleconetgermany
 twitter.com/medicaleconet

Distributor

*All specifications are subject to change without notice

ME-ENG-H25-ECTWLCD-01-17